



CareSens N POP Blood Glucose Meter

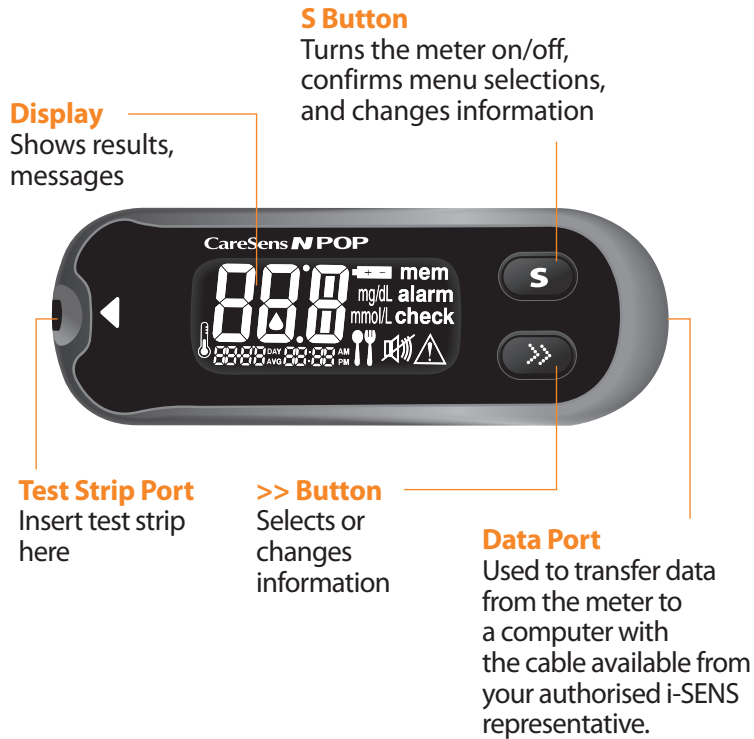

Note: The cable for data transmission to PC can be ordered separately. Please contact your authorised i-SENS sales representative.

CareSens N POP Blood Glucose Meter Display

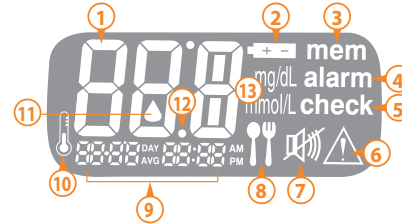

- 1 Test results :** test results displaying panel
- 2 Battery symbol :** indicates meter battery is running low and needs to be replaced
- 3 mem :** appears when test results stored in the memory are displayed
- 4 alarm :** appears when the post-meal alarm has been set
- 5 check :** appears when the control solution test results are saved or displayed
- 6 Hypoglycemia symbol :** appears when the test result is below the hypoglycemic level
- 7 Sound symbol :** appears only when the sound is set to OFF
- 8 Post-meal test flag :** appears during post-meal testing and when post-meal test results are displayed
- 9 Month / Day / Hour / Minute**
- 10 Temperature symbol :** displays recorded temperature when blood glucose levels are tested
- 11 Blood insertion symbol :** indicates meter is ready for the application of a drop of blood or control solution
- 12 Decimal point :** appears when the blood glucose measuring unit is set to mmol/L
- 13 mg/dL , mmol/L :** unit for measuring blood glucose

Note: The unit of measurement is fixed in mmol/L and it cannot be changed to mg/dL by the user.